

Certificate of Analysis

Specifications

Net peptide content	Nitrogen analysis
Purity Test	HPLC
Molecular Weigh	MS
Endotoxin	Endotoxin Analysis

Sample Summary


Product:	CP20	
Identity:	Confirmed	
Appearance:	White Lyophilized Powder	
LotNo.	N/A	
Reported:	2026.7.10	

Analytical Results

Test	Result
Identity (LC-MS)	N/A
Purity (HPLC-UV)	99.76%
Net Content	19.52mg
The endotoxin level in sample is	N/A

Method: HPLC with UV detection coupled with mass spectrometry (LC-MS).

Chromatogram

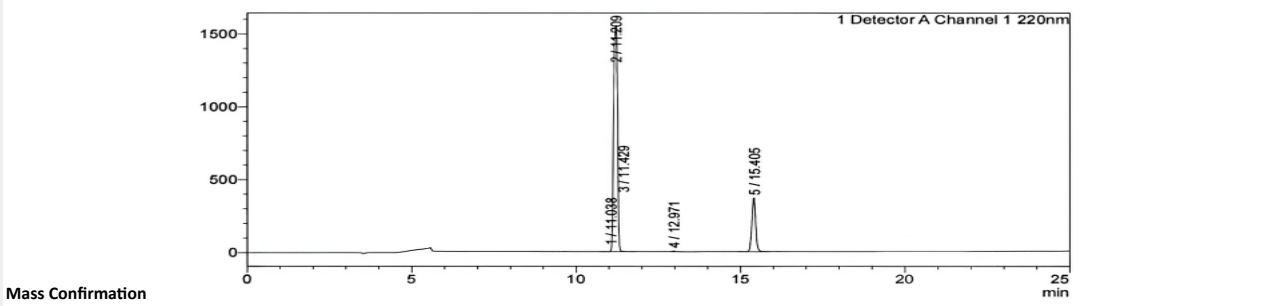

Mass Confirmation

Tested by	Miss.XieYuJing , QC	Date:2026/07/11
Tested by	Miss.ZhangBo , QC	Date:2026/07/11

The peptide purity analysis reported here was conducted using LCMS/MS under standard laboratory conditions. The net peptide content stated in this report is determined via nitrogen quantification using an elemental analyzer. This analysis is intended for informational purposes only and is specific to the sample(s) provided. The peptides tested are intended for research use only and are not approved for human or veterinary use, diagnostic, therapeutic, or clinical applications. Results should be interpreted by qualified professionals within the scope of the intended research. The accuracy and reliability of the test may be influenced by sample integrity, handling, and other experimental variables.

Test Conditions

HPLC

Pump A : 0.065% trifluoroacetic in 100% water (v/v)
Pump B : 0.05% trifluoroacetic in 100% acetonitrile (v/v)
Total Flow:1 ml/min
Wavelength:220 nm

<<LC Time Program>>

Time	Module	Command	Value
0.01	Pumps	B.Conc	5
25.00	Pumps	B.Conc	65
25.01	Pumps	B.Conc	95
27.00	Pumps	B.Conc	95
27.01	Pumps	B.Conc	5
35.00	Pumps	B.Conc	5
35.01	Controller	Stop	

<<Column Performance>>

<Detector A>

Column : Inertsil ODS-SP 4.6 x 250 mm
Equipment:ZJ19010014

MS

Interface	:ESI	Equipment	: ZJ17010408
Nebulizing Gas Flow	:1.5L/min	Interface Bias	: +4.5 kV
CDL Temp	:250°C	Drying Gas Flow	:5 L/min
Block Temp	:200°C	T.Flow	:0.2 ml/min
		B.conc	:50%H2O/50%MeOH

Endotoxin

LAL sensitivity: 0.02EU/ml

Sample Endotoxin level:

1. The endotoxin level in sample was detected using Horseshoe Crab Reagent, the phenomenon is listed as followed:

Table 1:

Inspection tube	0.1ml LAL+0.1ml LAL Reagent Water+0.1ml Detection solution	Parallel	
		1	2
Sample	0.1ml LAL+0.1ml LAL Reagent Water+0.1ml Sample solution	-	-
Sample positive control	0.1ml LAL+0.1ml LAL Reagent Water+0.1ml 2 λSample positive solution	+	+
Positive control	0.1ml LAL+0.1ml LAL Reagent Water+0.1ml positive solution	+	+
Negative control	0.1ml LAL+0.1ml LAL Reagent Water+0.1ml LAL Reagent Water	-	-

Nitrogen analysis

Calculation Method:

The net peptide content/%= $N\% \times MW(\text{peptide}) / (n \times MW(\text{nitrogen}))$

N%: Nitrogen content determined by elemental analysis;

MW(peptide): Theoretical Molecular Weight;

n: The number of nitrogen element;

MW (nitrogen): atomic weight(nitrogen);